

HAIYE




HP Series Reusable Anesthesia Fiber Optic Laryngoscope HP-W/HP-Z

Application

It is used clinically to provoke the patient's epiglottis to expose the glottis, and to guide medical staff to perform tracheal intubation and mechanical ventilation.

Suitable for anesthesia departments, ICU, pediatrics, emergency centers, ambulances, etc. in medical institutions.

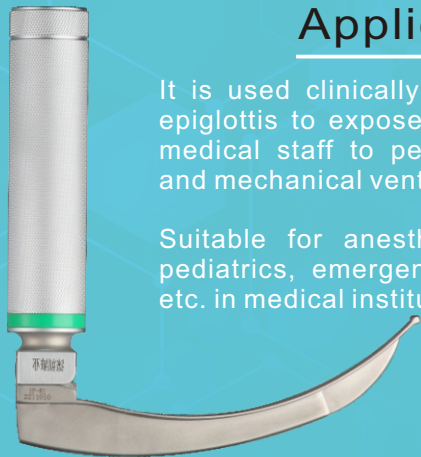

Feature

This product is divided into two series, HP-W and HP-Z, corresponding to adults and children respectively. Using high brightness LED cold light bulb, 4mm fiber bundle, high brightness, wide vision, built-in fiber optic design to avoid breakage and wear, convenient for cleaning and disinfection, and extend the service life; The laryngoscope blade is made of stainless steel and cast in one piece, which is made in one-step forming, using a three-point fixing method.



Parameter

Model of laryngoscope blade	Composition	Size of laryngoscope blade	Applicable user
HP-W1	Mac 4	155*15 mm	Adult
HP-W2	Mac 3	134*14.2 mm	Adult
HP-W3	Mac 2	114*14 mm	Adult
HP-Z1	Mac 1	93*10.9 mm	Pediatric
HP-Z2	Miller 1	80*9.4 mm	Pediatric
HP-Z3	Miller 0	68*10.7 mm	Pediatric



Shenzhen Haiye Medical Technology Co., Ltd.
Add: Grand Equipment Industrial Park, Shenzhen, Guangdong, China
Tel: +86-29-87872952
Email: info@haiyemedical.com
Website: www.hyendoscope.com