

HAIYE



## HG5.2 Video Stylet



### Application

It is provided for clinical use to provoke the patient's epiglottis to expose the glottis, and to guide medical staff to perform tracheal intubation and implement mechanical ventilation.


Suitable for anesthesia departments, ICU, pediatrics, emergency centers, ambulances, etc. in medical institutions



### Features and Configuration

Product structure: handle and insertion part

This product can be used with 3.5" touch screen display or 12.1" (non-touch), 13.3" big touch screen monitor.

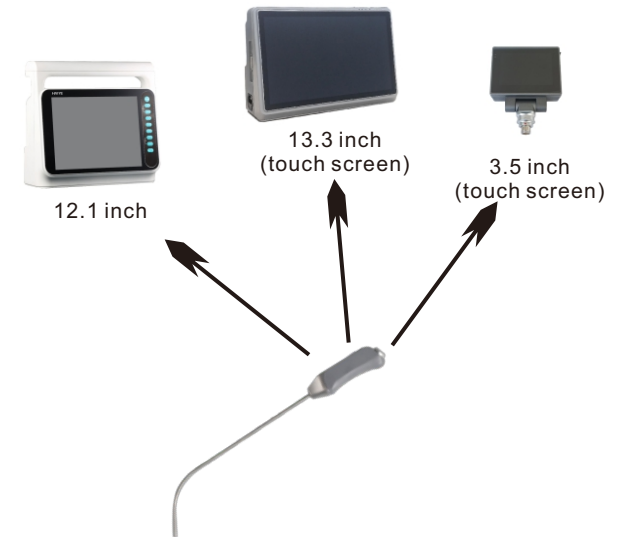

Field angle	70°±9°
Feild depth	3~100mm
Size	diameter 5.2±0.1mm, length 400±10mm, bending angle 110°±10°
Camera	2 million pixels; resolution ≥7.13lp/mm
Illumination	≥3000lx
Function	connect to 3.5", 12.1" or 13.3" display, one-click photo and video



**Shenzhen Haiye Medical Technology Co.,Ltd.**

Add: Grand Equipment Industrial Park, Shenzhen, Guangdong, China

Tel: +86-29-87872952

Email: info@haiyemedical.com

Website: www.hyendoscope.com