




HP-Y

Single-use Anesthesia Fiber Optic Laryngoscope

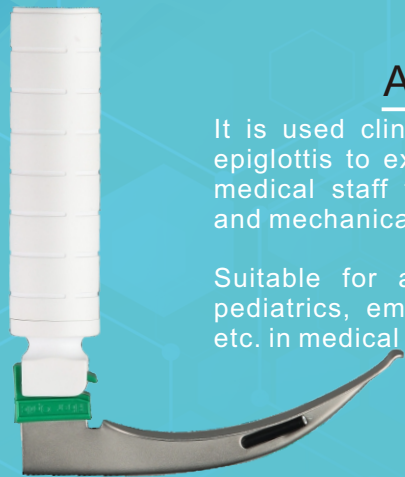

Application


It is used clinically to provoke the patient's epiglottis to expose the glottis, and to guide medical staff to perform tracheal intubation and mechanical ventilation.

Suitable for anesthesia departments, ICU, pediatrics, emergency centers, ambulances, etc. in medical institutions.



Feature and Configuration

The entire product is packaged in sterile packaging and is for one-time use, requiring no cleaning or disinfection; The laryngoscope blade is made of stainless steel and is available in 6 specifications, which can meet the needs of different groups of people; The handle has a built-in high-brightness LED cold light source bulb, which is transmitted to the front end through the optical fiber bundle on the laryngoscope blade. High brightness, wide field of view.

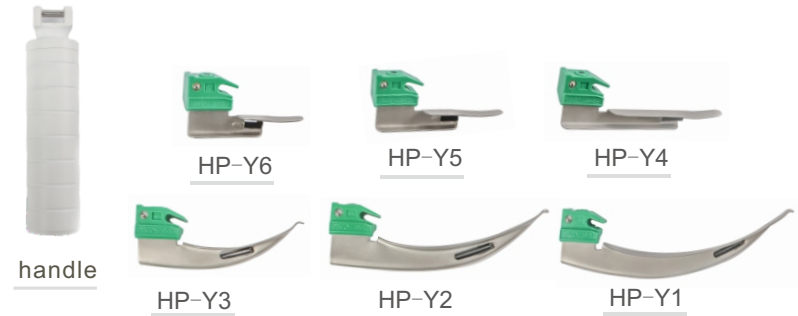

Parameter

Model of laryngoscope blade	Shape of laryngoscope blade	Size of laryngoscope blade	Applicable user
Y1	Curved	153*13.8 mm	Adult
Y2	Curved	132*14.8 mm	Adult
Y3	Curved	100*13.8 mm	Adult
Y4	Straight	101.7*9.7 mm	Pediatric
Y5	Straight	80*10.2 mm	Pediatric
Y6	Straight	67*10.5 mm	Pediatric



Shenzhen Haiye Medical Technology Co., Ltd.
 Add: Grand Equipment Industrial Park, Shenzhen, Guangdong, China
 Tel: +86-29-87872952
 Email: info@haiyemedical.com
 Website: www.hyendoscope.com